


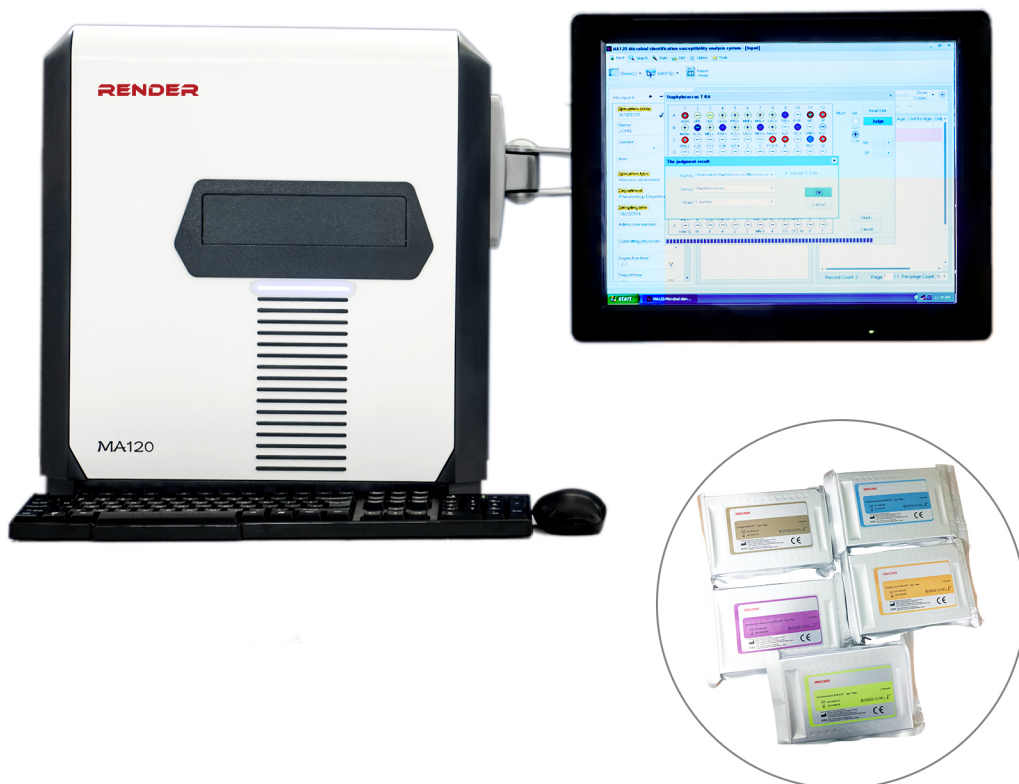


MA120

Automated ID&AST System

- A system comprised of a PC & Software and Reader, affordable for small & medium size laboratories to do automatically identification and antibiotic susceptibility tests.
- Capable of performing Identification and Antibiotic Susceptibility tests for Nonfermenter, Streptococcus/Enterococcus, Staphylococcus, Fungus and Enterobacteriaceae. It covers more than 300 species bacteria and fungi.
- The Test Panel consists of 120 microwells, which integrates 24 microwells for Identification and 96 microwells for Susceptibility test. Both results can come out in a single test.

● Specifications

Capacity: less than 1 minute/test panel
 Data Interface: USB, Bidirectional interface
 available for LIS or HIS
 Input: 100-240 VAC, 50/60HZ
 Dimension: 400*390*440mm
 Weight: 14.5kg

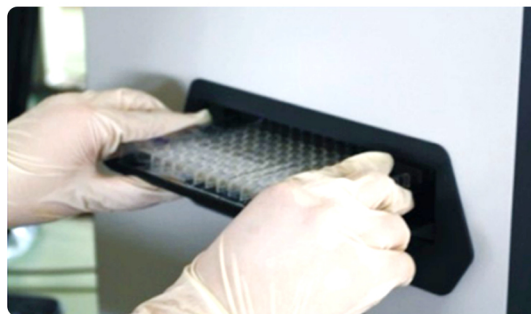


● Identification

- Detection Principle: Colorimetric Method
- Visual inspection is applicable


● Susceptibility Test

- Detection Principle: Turbidimetric Method
- True MIC for Susceptibility test, Unique Patent design microwells
- Up to date antibiotic susceptibility test and grouping per CLSI/EUCAST standard

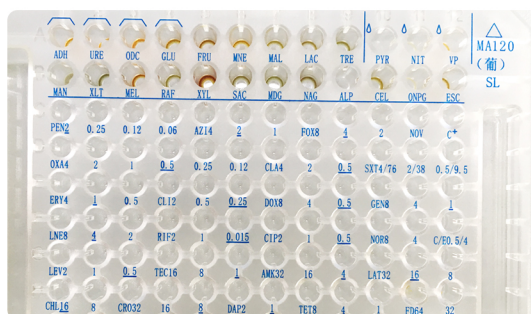

● Software

- User-oriented design
- Expert System included, marking antibiotic choice easy by glance
- Independent quality control for each test panel
- Upgrade annually according to CLSI/EUCAST
- Hospital infection system, case management system and statistical system make your work simple and efficient

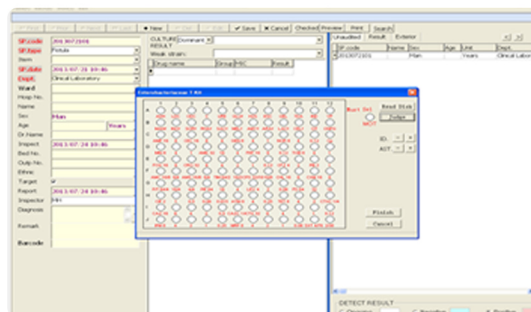

● Flexible and More Options for Automation



Turbidimeter(MT-06) and Automatic Sampler(AS120) are optional to ID&AST system.

Render Biotech Co., Ltd.

Add: 6K#, Building No.1, Baisha Emerging Industry Zone, Shahexi Road, Nanshan District, 518052, Shenzhen, China

Tel.: +86 755 8653 6648 E-mail: info@render-health.com

Web: www.renderbio.com